

HERO NXT

The "hero" for the little ones

Hernik



~ 3 months -
4 years

61 - 105 cm

up to 18 kg



Comfortable & secure

We rely on the latest technology to make our products safe.

- ✓ Best ratings on comparison portals
- ✓ Approved acc. to R129



Trust in experience

20 years of experience developing and manufacturing the best products.

- ✓ The only rear-facing car seat for small children with full European approval



Made in Germany

All our products are developed and manufactured in Germany.

- ✓ Certified according to ISO 9001

Maximum protection for the little ones

The HERO-NXT is a car seat modified for the rehabilitation market. Thanks to the 360° rotation function, the seat can be used facing backwards and forwards. Due to the weak neck muscles of children under 15 months, the HERO-NXT should be used in a rearfacing position.



1 5-point harness

The integrated 5-point harness secures your child in the car seat. The quick-release fastener makes it easy and quick to put your child in and take him out.

2 Inner upholstery

The inner upholstery developed by HERNIK ensures optimal positioning and significantly more comfort, especially for small children

3 Protective bracket and carrying handle

The carrying handle helps you to carry the child seat ergonomically and safely. While driving, it also acts as a protective bracket and ensures stability in the event of a rear-end collision.


4	Total height	66 - 73 cm
5	Total width	50 cm
6	Total depth	65 cm



Quick installation

The child car seat is quickly and easily fixed in the car using the ISOFIX connection. It is not necessary to secure the child with the 3-point belt.

An ISOFIX system is a prerequisite for our HERO-NXT, which is why the connection is already integrated in the basic equipment.

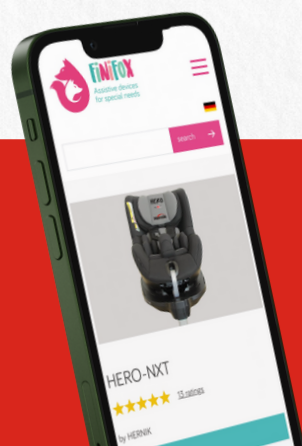

We are very proud of the many positive customer ratings. On our website you can find out, why HERNIK regularly achieves top ratings on comparison portals.

Markus Jacksch,
Managing Director



Individual accessories for every child

1 2 Adjustable back height/ Adjustable seat angle

To support children with limited upper body control, the seat angle can be adjusted, allowing you to switch between an active and a more passive position. In addition, the child car seat grows with your child thanks to the adjustable headrest.

3 Side supports

Children with poor trunk control usually need a firm, stable hold. This is made possible by the aluminum side supports, which can be adjusted to the child's body in width and height with complete flexibility. This absorbs unavoidable forces when accelerating, braking and cornering.

4 Security tray

The security tray provides more stability and a more upright posture when sitting.

5 Abduction block

To prevent dislocation of the hips and compression of the knees (due to muscle tension), an abduction block can be chosen.

6 360° rotatable

The HERO-NXT can be rotated 360° and can be used both in and against the direction of travel. For easier transfers, the child car seat can be swiveled 90° so that you can put the child in the car seat without overstraining your back. Buckle-in is also much easier in this position. Then you rotate the HERO-NXT into the desired position.

7 Adjustable support leg

Using the adjustable support leg, the HERO-NXT can be easily installed in the vehicle. The support leg gives the child car seat additional stability above the vehicle floor and at the same time minimizes the tipping movements of the seat.

We are certain that we can only develop the best seats together with our customers. That is why we are very happy to receive your suggestions for improvement and will incorporate them into our development.

Please send us a message:
office@hernik.de



Further products from HERNIK



Child car seats for every size

HERNIK accompanies your child from infancy to adulthood (from 0-75 kg). For more complex requirements and individual requests we also offer customised solutions.

Head protection helmets and seat shells

With our wide range of head protection helmets we offer every user the necessary protection against injuries. Our CAPOVARIO in particular is a real highlight that makes everyday life much easier thanks to its special design.

We also supply you with individual seat shells, such as the ANATOMIC UNICO and ANATOMIC FLEX aluminium shells.



Hernik

**Discover all our products
on hernik.de**

HERNIK GmbH

Gewerbering 8,
83543 Rott am Inn

+49 (08039) 9093930
office@hernik.de
www.hernik.de

