

IPAI LGT

light, safe and comfortable



 ~3.5 - 12 years

100 - 150 cm

up to 36 kg



Comfortable & secure

We rely on the **latest technology** to make our products safe.

- ✓ **Best ratings** on comparison portals
- ✓ **Approved acc. to R129**



Trust in experience

20 years of experience developing and manufacturing the best products.

- ✓ **Rehab car seat with full European approval**



Made in Germany

All our products are developed and **manufactured in Germany.**

- ✓ **Certified according to ISO 9001**



The lightweight car seat

The IPAI-LGT is a modified car seat from Britax Römer, which is characterized by its low weight and an extensive basic equipment.


The seat comes standard with a 5-point harness, a large upholstery for the belt buckle, an ISOFIX adapter and side supports. The standard equipment also includes SICT side impact protection. This increases safety in the event of a side impact without you having to sacrifice space on the back seat.

Security features

- 1 SICT side impact protection
- 2 Side supports
- 3 Shoulder pads
- 4 5-point harness
- 5 Upholstery for the belt buckle
- 6 Headrest multi-height adjustable

Dimensions

- | | | |
|---|--------------|------------|
| 7 | Total height | 67 - 85 cm |
| 8 | Total width | 48 cm |
| 9 | Total depth | 39 cm |

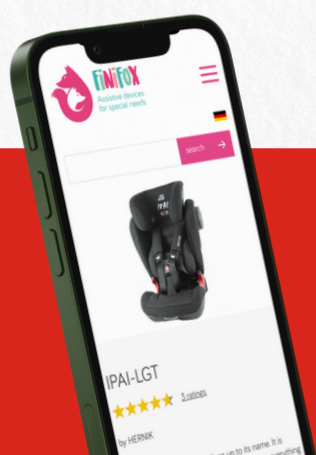

We are very proud of the many positive customer ratings. On our website you can find out, why Hernik regularly achieves top ratings on comparison portals.

Markus Jacksch,
Managing Director

As individual as your child

Innovative accessories allow you to adapt the seat to your child's individual needs. It is great, that all accessories can be retrofitted in just a few simple steps. These include:

① **Seat depth adapter**

The seat depth adapter allows you to extend the seat depth from +2 to +12 cm. This improves comfort and allows the seat to grow with your child.

② **Seat angle adjustment**

The seat angle adjustment allows for a more comfortable sitting position, which increases well-being, especially on longer journeys. This seat inclination may also be useful for children with limited body tension.

③ **Security tray**

The security tray provides more stability and a more upright posture when sitting.

④ **Abduction block**

To prevent dislocation of the hips and compression of the knees (due to muscle tension), an abduction block can be chosen.



Easy to use via ISOFIX

The ISOFIX system in cars consists of metal hooks on the vehicle seat and brackets on the child seat. When the child seat clicks into place, a secure connection is created that increases stability and safety, especially in the event of an accident.

Using a central handle on the front, the depth of the ISOFIX adapter can be adjusted with one hand. This allows the seat to be pulled forward if necessary, which opens the back angle and allows a resting position to be quickly adopted. If the seat is used in different vehicles, the ISOFIX system allows you to install the seat quickly and safely. If the seat is not needed, the ISOFIX system disappears into the seat unit.



We are certain that we can only develop the best seats together with our customers. That is why we are very happy to receive your suggestions for improvement and will incorporate them into our development.

Please send us a message:
office@hernik.de



Further products from HERNIK



Child car seats for every size

HERNIK accompanies your child from infancy to adulthood (up to 75 kg). For more complex requirements and individual requests we also offer customised solutions.

Head protection helmets and seat shells

With our wide range of head protection helmets we offer every user the necessary protection against injuries. Our CAPOVARIO in particular is a real highlight that makes everyday life much easier thanks to its special design.

We also supply you with individual seat shells, such as the ANATOMIC UNICO and ANATOMIC FLEX aluminium shells.



Hernik

**Discover all our products
on hernik.de**

HERNIK GmbH

Gewerbering 8,
83543 Rott am Inn

+49 (08039) 9093930
office@hernik.de
www.hernik.de

