

KIDSFLEX^{PRO}

Hernik

Clever design - for what matters most!



Pro 1: ages 3 to 9 years

Pro 2: ages 6 to 14 years



Comfortable & secure

We rely on the **latest technology** to make our products safe.

✓ **Fast switch** between active and passive **seating position**



Trust in experience

20 years of experience developing and manufacturing the best products.

✓ Rehab car seat with **full European approval**



Made in Germany

All our products are developed and **manufactured in Germany.**

✓ **Certified according to ISO 9001**

Simply outstanding

The KIDSFLEX-PRO is a seat engineered for rehabilitation, allowing for optimal adaption to a wide range of individual needs. Even in its standard configuration, it comes equipped with a 5-point harness, a **width-adjustable headrest (1)**, thorax pads, shoulder supports, a seat-depth extension, back-angle adjustment and an ISOFIX adapter.

The newly designed **rear section (2)** can be adjusted telescopically, enabling simultaneous changes to the positions of the headrest, shoulder straps and thorax pads.

In addition, the KIDSFLEX-PRO is equipped with a **seat-depth extension (3)** that increases the seat depth by up to 6 cm.



Rapidly adaptable

A completely redesigned **ISOFIX quick-adjustment system (4)** allows you to switch between active and passive seating positions within seconds. It creates the necessary space between the KIDSFLEX-PRO and the vehicle seat backrest to adjust both the **backrest angle (5)** and the **seat angle (6)**. And all of this is done quickly and easily — without any tools.

The new KIDSFLEX-PRO has successfully passed extensive crash tests in accordance with the R129 standard and has obtained full European approval with an E-number.



Simple and effortless adjustments —
no tools necessary!

We are very proud of the many positive customer ratings. On our website you can find out, why Hernik regularly achieves top ratings on comparison portals.

Markus Jacksch
CEO

PRO-Highlights

With the **PRO seat cushion (1)**, we introduce an ergonomically designed cushion that provides optimal support and significantly enhances seating comfort. The integrated, width-adjustable leg guides help maintain a stable leg position, while the abduction padding promotes a healthy, ergonomic sitting posture. This ensures that your child receives the best possible support, contributing to greater safety, comfort, and well-being on every journey.

The **PRO headrest (2)** offers increased safety compared with the standard headrest. It features side impact plates and additional adjustment options enabled by patented predetermined bending points.

The newly developed **rotating mechanism (3)** provides essential support for parents. The ability to rotate the seat by 90° greatly simplifies placing the child into the seat and securing them.



Accessories

For optimal positioning, additional accessories such as various **thorax-pads (4)**, a **footrest (5)**, an **abduction wedge (6)** and a **security tray (7)** are available.

Please refer to our order form for a complete overview of all available options.



We are certain that we can only develop the best seats together with our customers. That is why we are very happy to receive your suggestions for improvement and will incorporate them into our development.

Please send us a message:
office@hernik.de



Further products from HERNIK




Child car seats for every size

HERNIK accompanies your child from infancy to adulthood (up to 100 kg). For more complex requirements and individual requests we also offer customised solutions.

Head protection helmets and seat shells

With our wide range of head protection helmets we offer every user the necessary protection against injuries. Our CAPOVARIO in particular is a real highlight that makes everyday life much easier thanks to its special design.

We also supply you with individual seat shells, such as the ANATOMIC UNICO and ANATOMIC FLEX aluminium shells.

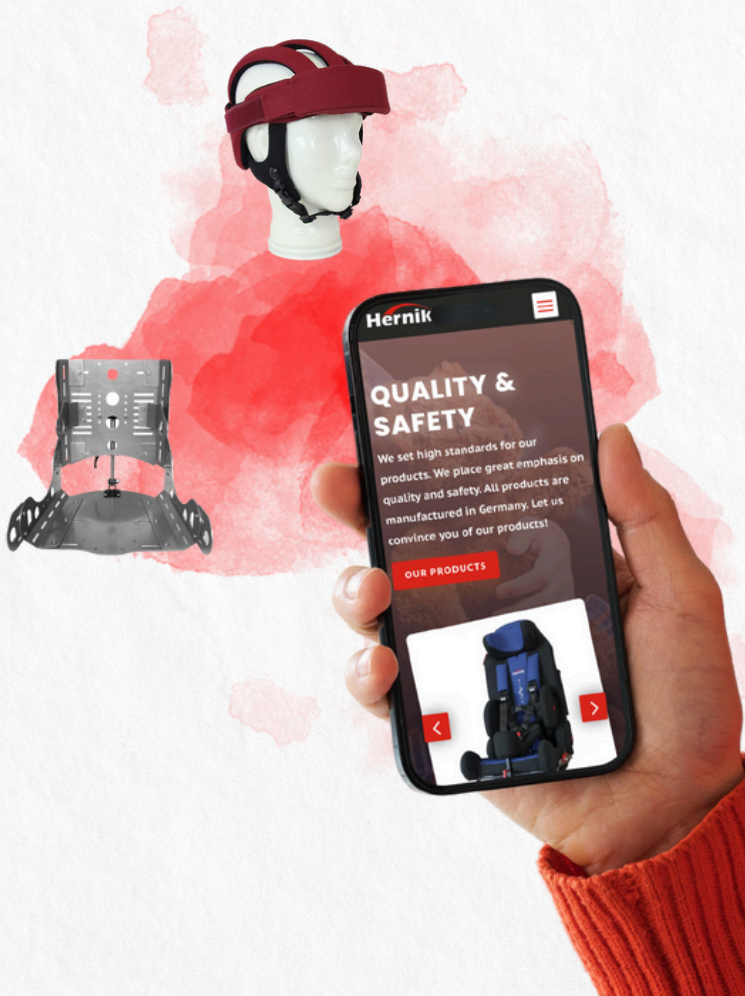

Hernik

**Discover all our products
on hernik.de**

HERNIK GmbH

Gewerbering 8
83543 Rott am Inn

+49 (08039) 9093930
office@hernik.de
www.hernik.de

